


Babytherm 8000 OC with Radiant Heater RH600 Technical Service Bulletin # 8

Re.: Replacement of swivel bearing 2M 18 652 of RH600

Update: August 24, 2000

Reference Doc: N/A

Reason: To ensure the swivel bearing version with one nut does not become loose.

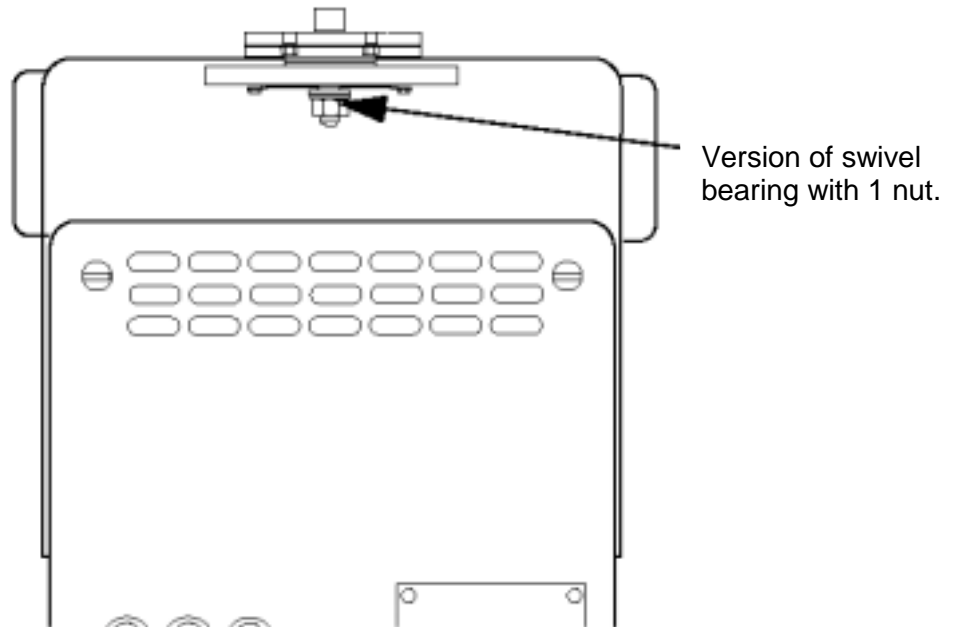

Solution: Installation of a newer version with 2 nuts P/N 2M 18 652 or P/N SE 2M 18 652.

Devices affected: Radiant Heaters RH600 in the USA and Canada shipped prior to April 1996 up to S/N ARKD-0040 (not all devices are affected).

Please note: The RH600 could be delivered as a stand-alone unit and with the following combinations:

RH600 with wall mount	P/N 2M 20 116
RH600 with mobile stand	P/N 2M 20 539

Babytherm 8000 OC

P/Nos. FR00048, FR00050 or
FR00086

When: At next Service Call or Preventive Maintenance in line with the actions of TSB Radiant Heater RH600 # 5.

Additional Info: Modified units RH600 are labeled with "New Sensor Technology" on the keypad.

[This TSB is identical to TSB Radiant Heater RH600 # 4.](#)

Cost: Parts free of charge.

Ordering Info: 2M18652 or SE 2M18652

Distribution: Dräger Service Personnel, Authorized Service Organizations, and trained hospital technicians for Radiant Heater RH600.

If you have any questions, please contact Technical Support by phone at 1-800-543-5047 or by fax at 1-215-723-5935

Dräger Medical, Inc.
Technical Product Manager